

LIMITED WARRANTY

The Company will repair this product with new or rebuilt parts, free of charge for one year from the date of original purchase to the event of defect in material or workmanship.

A purchase invoice or any other proof of date of original purchase will be required before warranty performance is rendered. This warranty does not cover any damage which results from accident, misuse, abuse, alteration, modification or commercial use of this product.

This warranty is not available if the product is used in a health club or for any commercial purposes.

CUSTOMER NAME: _____

DATE: ___/___/___

DEALER NAME: _____


DEALER ADDRESS: _____

tokuyo

INSTRUCTION MANUAL

CONVERTIBLE BEAUTY LEG MASSAGER

MODEL NO.: TF-655

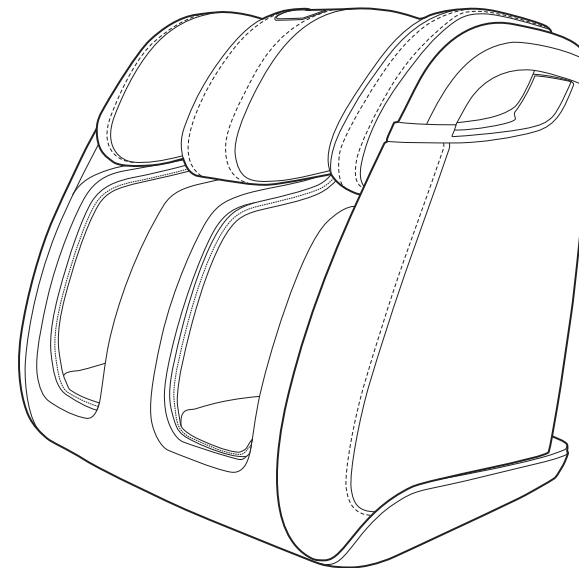

Thank you for purchasing our Convertible Beauty Leg Massage.
Be sure to read this instruction manual carefully and follow it properly.
Keep this manual within reach at all times for future reference.

Thank you for purchasing our Convertible Beauty Leg Massager. Be sure to read this instruction manual carefully and follow it properly. Keep this manual within reach at all times for future reference.

CONTENTS

IMPORTANT SAFETY INSTRUCTIONS.....	02
NAME OF EACH PARTS.....	06
FUNCTION AND CONTROLS.....	08
BEFORE USEAGE.....	09
REMOTE CONTROL INSTRUCTIONS	10
AFTER USEAGE	12
CARE AND MAINTENANCE	13
Q&A	14
SPECIFICATION	14

TROUBLE SHOOTING

Symptoms	Possible causes	Solution
Product not running	Power cord not properly plugged.	Check if product's Power Supply Voltage is compatible with your country's power supply. And make sure power cord is plugged correctly.
	15 minutes over. automatically switch to standby mode.	If you wish to have massage one more time, press ON/OFF button again.
	The leg massager automatically stops after long period of use more than 30 min to prevent overheating.	Before restart the product, must be cut off the power supply and let it cool down for 60 minutes.
Overheating	The leg massager auto stops after long period of use more than 30 min	Cut off the power supply and let it cool down.
Makes sound when running	Airbags on product will cause sound while working.	This sound is normal, does not indicate problem with product.

SPECIFICATION

Model No.	TF-655
Rated Voltage and Frequency	<input type="checkbox"/> 100~120V <input type="checkbox"/> 220~240V <input type="checkbox"/> 50Hz <input type="checkbox"/> 60Hz
Rated Power	80W
N.W. / G.W.	17.2 / 19.8 kg
Dimension	510(L)*455(W)*497(H)(mm)
Origin	China
Serial No.	


MAINTENANCE




WARNING

UNPLUG LEG MASSAGER BEFORE CLEANING.
DO NOT TOUCH POWER CORD WITH WET
HANDS.


Wet a soft towel in low concentration kitchen solution, wring out thoroughly. Wipe leg massager with wet towel.

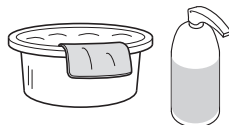

Use a dry towel to wipe leg massager dry afterwards.

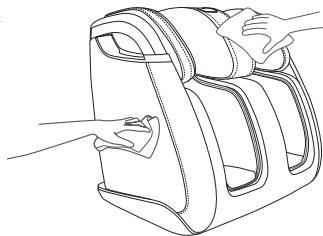

Lift up the Convertible leg Part.
Find zipper slider of cover
Unzip cover and take out fabric cover
Apply cover back on after cleaned and dried.

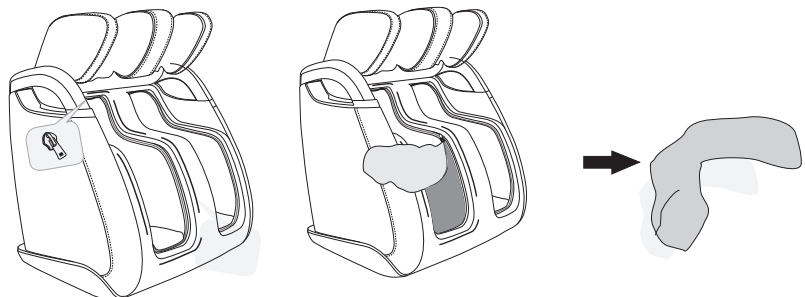

Keep away from cleaning dissolvant Will endanger product to color fading, deformation or other damages.
Remove fabric cover carefully.
If the fabric is torn, do not use leg massager, contact customer service immediately.
Using Leg massager without the fabric cover will cause serious danger.
Please avoid at all costs.
Do not let children play with or near Leg Massager while cover is off.


IMPORTANT SAFETY INSTRUCTIONS

Thank you for purchasing Convertible Beauty Leg Massager!
Before using this product, read IMPORTANT SAFETY INSTRUCTIONS carefully to enjoy leg massager safely and prevent accidents.

IF YOU MEET THE FOLLOWING CONDITIONS, PLEASE CONSULT A DOCTOR BEFORE USE.

- *If you are sensitive to electromagnetic wave; or have an implanted medical device.
- *If you suffer from a malignant tumor.
- *If you suffer from heart disease.
- *If you suffer lack of temperature sensitivity.
- *During menstruation period or during pregnancy.
- *If you suffer from diabetic neuropathy.
- *If you have an injury.
- *If you are highly sensitive to noise.
- *If you are having a fever.
- *If you are feeling weak.

- *Use correct input voltage (Not following this instruction might cause fire or electricity shock).
- *Check for rips and tears on the fabric cover. If damage is found, no matter how small, unplug power cord immediately and send for repair.
- *When moving leg massager make sure your surroundings is free of people, pets or other obstacles.
- *After using leg massager, switch power off, this will lock the rotation systems and prevent product from hurting pets and children.
- *Do not sit on a wheeled chair while using leg massager.

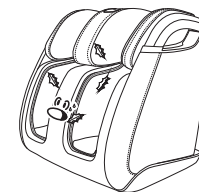

IMPORTANT SAFETY INSTRUCTIONS

*Do not use for a long period of time, might cause cryogenic burns

*This product is design for leg massaging only. Do not use to massage other body parts.

*Do not let children use leg massager. When leg massager is running keep children away.


*Do not stand or sit on leg massager

*Do not step on product forcefully

*Keep hands away from massaging mechanism to prevent getting stuck

*Non professional technicians are not allowed to dismantle or repair product.

*Keep power cord unobstructed at all times.

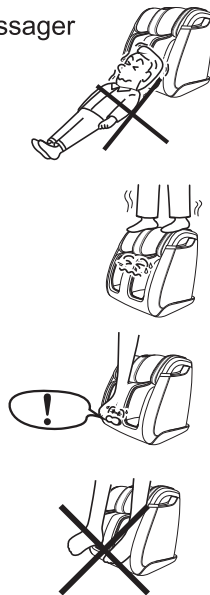

AFTER USAGE

After use



WARNING


Turn the power switch OFF and unplug the power cord after use.



ATTENTION

Unplug the power plug when not using.


Set leg massager to default position

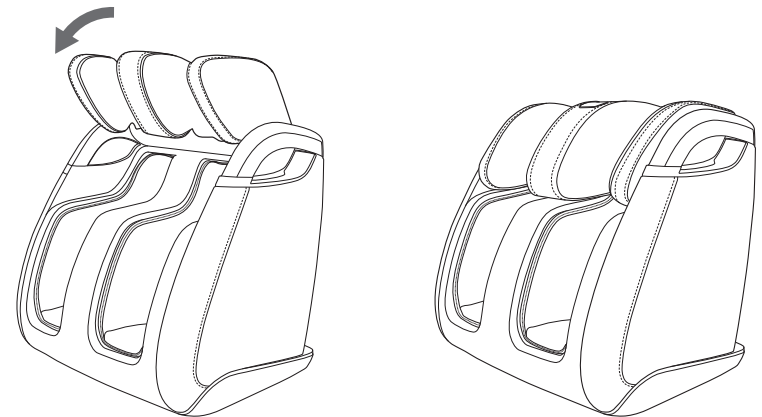

MANUAL MODE

There are 3 types of manual modes.
Press manual mode button to change to another mode.

Light of current working mode will be shining on the panel.

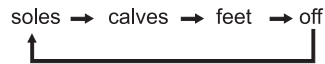

AIRBAG STRENGTH ADJUSTING BUTTON

6 strength levels to choose, default set is medium strength

Strength level is shown by the light signals in the manner seen below:



LO MID HI



HEAT BUTTON

Press heat button to turn on heat mode, press again to turn off



GENTLE MODE BUTTON

When foot massage feature is on, press this button to set bottom massaging strength lower.



KNEADING MODE BUTTON

Press to turn on kneading mode, press again to turn off



IMPORTANT SAFETY INSTRUCTIONS

*Do not wrap the power cord around the leg massager.

*Non-professional technicians are not allowed to dismantle or repair product

*When cleaning product, make sure power cord is unplugged.

*Do not touch power cord with wet hands.

If you feel sick during usage, stop process immediately and consult doctor


*Do not use product for a long period of time, take a timeout between sessions.

If leg massager does not startup, to avoid accidents, unplug power cord and contact your local supplier for inspection and repair.

*When unplugging do not pull by the cord. Clean power plug periodically to avoid fire due to dust.

*If leg massager was left unused for a long time, before usage, make sure it runs properly and the power cord is in good condition.

*Do not store product in damp places. Exposing product to dust and humid environment for long time period might cause danger and damage.

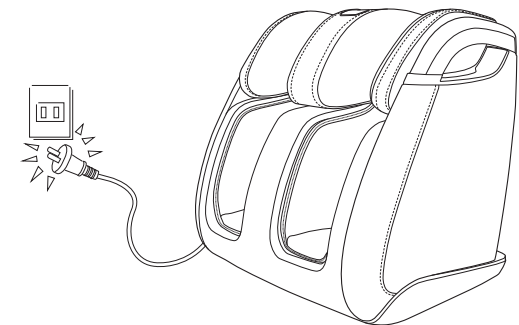

IMPORTANT SAFETY INSTRUCTIONS

*Do not use this product if your feet are smaller than 21cm.

*Do not use product in bathrooms or other humid places.

If an electricity blackout occurs, unplug power cord, this will prevent unwanted accidents when electricity is back on.

Connect to power, Switch ON/OFF to ON.

When power light is on, press power button on the control panel.



Power button

Emergency pause button

Pressing this button during a massage session will stop the process



Pause button

Auto mode

This product comes with 3 different auto massage modes

All modes are designed according to scientific studies.

Press Auto button

There are 3 types of auto modes, SHIATSU RELAX and QUICK.

Press auto button to switch mode.



- SHIATSU
- RELAX
- QUICK

Auto mode





Before use

*Product working environment
Place leg massager on a leveled surface


*Keep at least 1m away from televisions radios and other av appliances.
Doing this will avoid disturbing noise.

*Check fabric cover.
Before use, make sure fabric cover is free of rips and tears.
If damage is found, no matter how small, unplug and send for repair.
Power cord must be securely plugged


*Put your foot in leg massager
Place heel at the back of "legs" part of leg massager
adjust leg massager to preferred angle

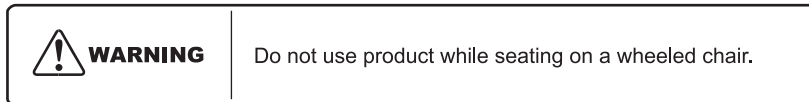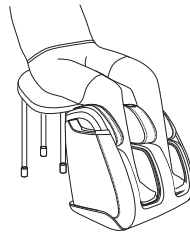

Default leg angle

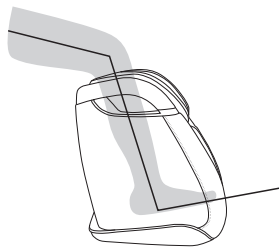


Lift up Convertible Leg Part to give knees a hot thermal patch treatment

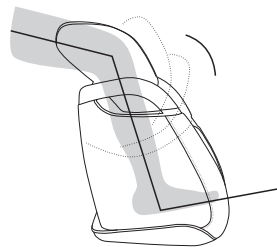

Adjust by hand the desired angle of Convertible Leg Part

Names of Functions

Front

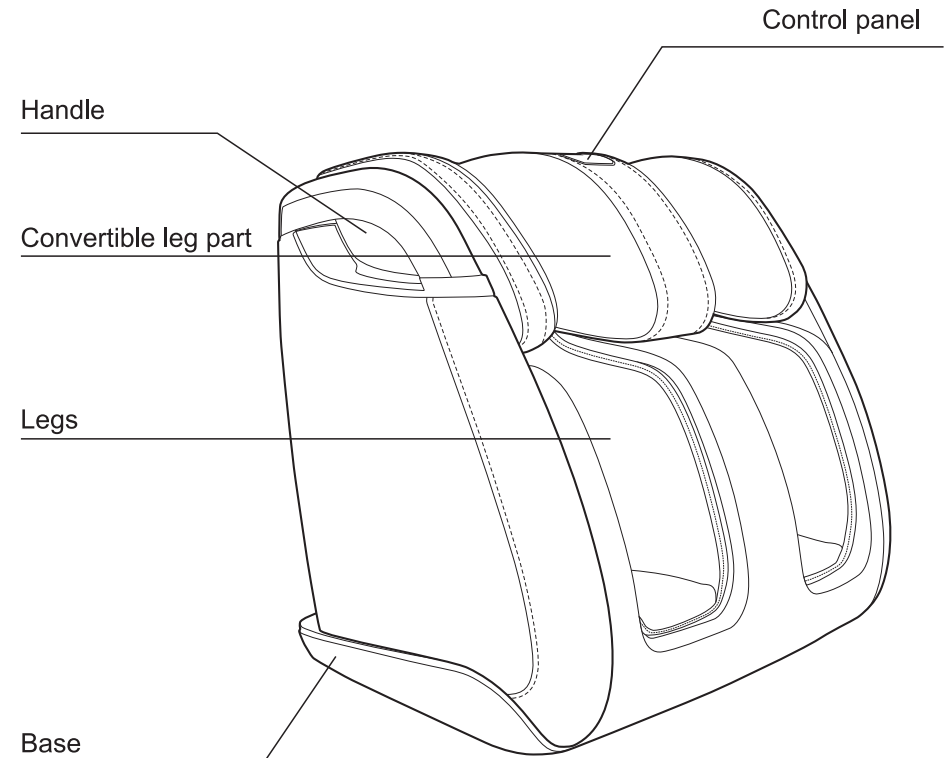

Names of Functions

Remote Control

Back

